


Quick start guide single-cell analysis ImageJ plugin

Installation of ImageJ and plugin.....	2
Image analysis of raw data	4
Pre-processing	4
Post-processing	7
Measurement data.....	11
Exemplary data sets	11
Troubleshooting	11


Installation of ImageJ and plugin

- Download current version of FIJI
<http://jenkins.imagej.net/job/Stable-Fiji/lastSuccessfulBuild/artifact/fiji-win64.zip>
- Unzip file
- Lunch fiji-win64.exe
- Add repository of the plugin to the ImageJ Updater as following:

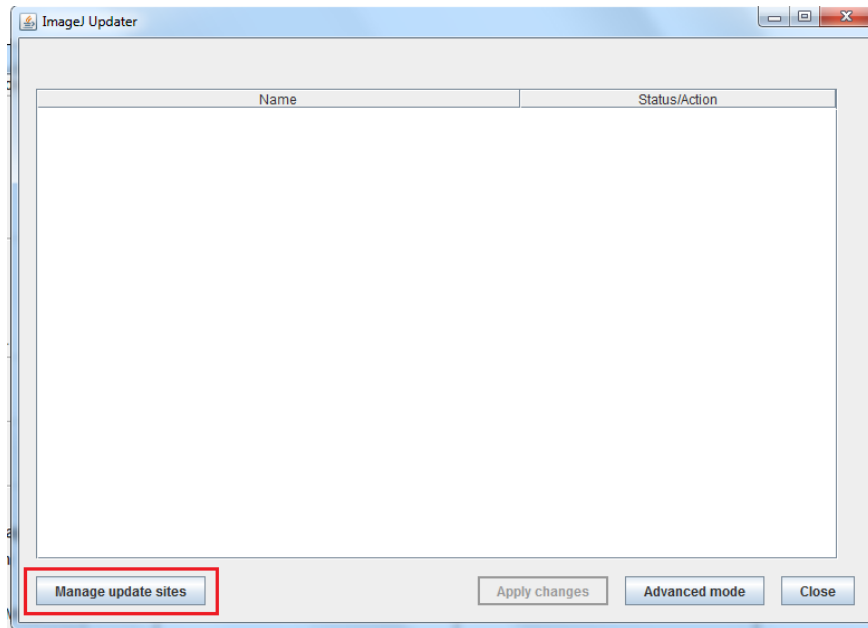

ImageJ Updater window

- Go to *Help->Update Fiji* and wait till the ImageJ Updater windows opens.
- Press the button "*Manager update sites*"

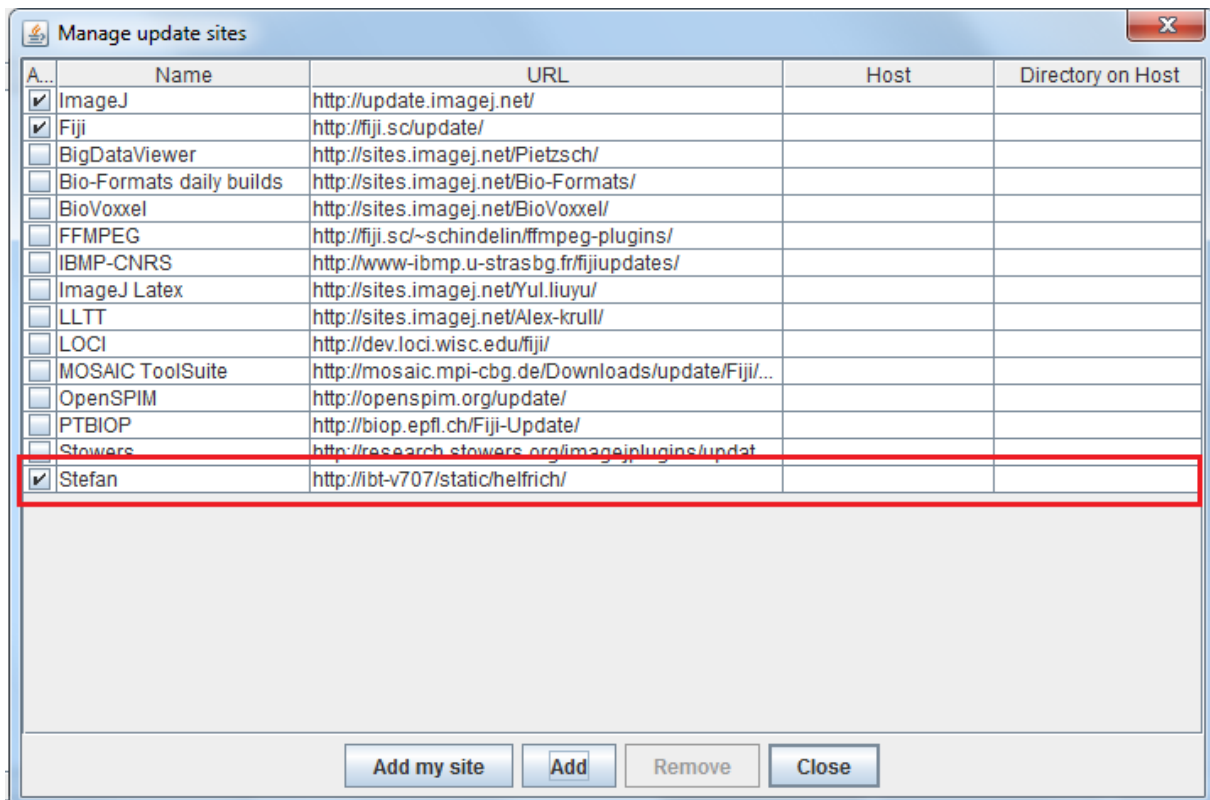


- Add the URL “http://ibt-v707/static/helfrich/” to the update sites using *Add*
- Press Close afterwards and Fiji will install the plugin


Image analysis of raw data

Pre-processing

- Export selected ND file as tif image series

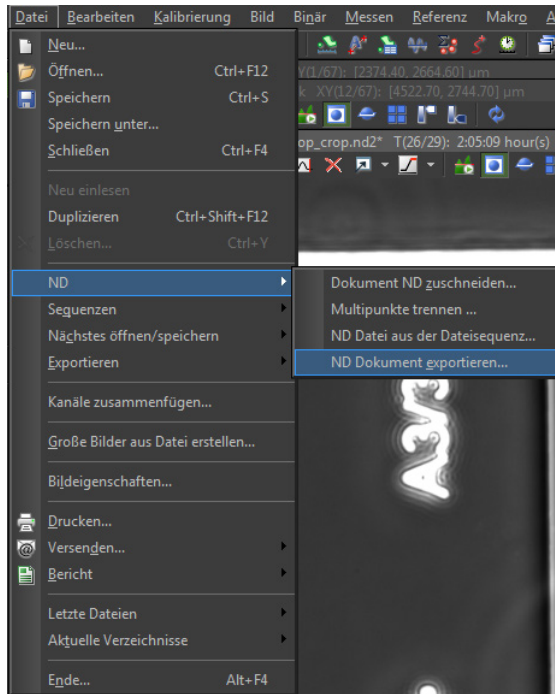


Export ND times series in NIS Elements

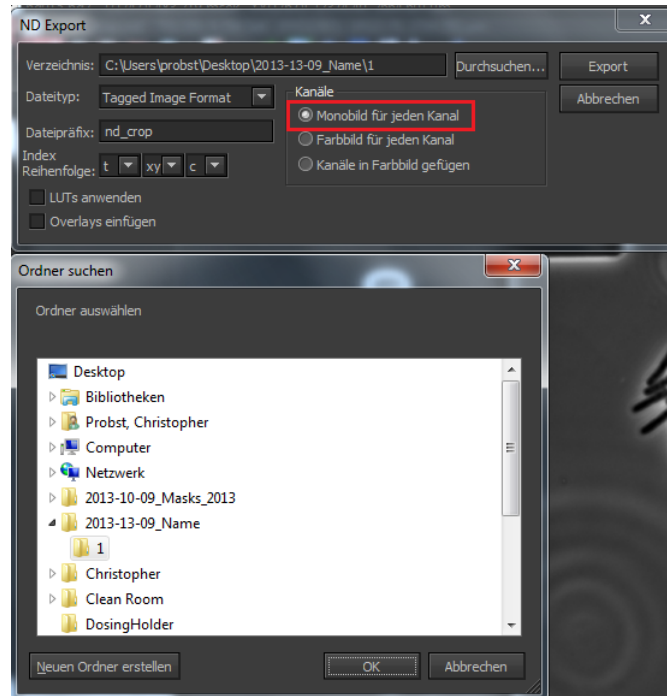


Select a folder for exporting the tif images


- Lunch Fiji if not done so already
- Import tif image series *File->Import->Image Sequence*

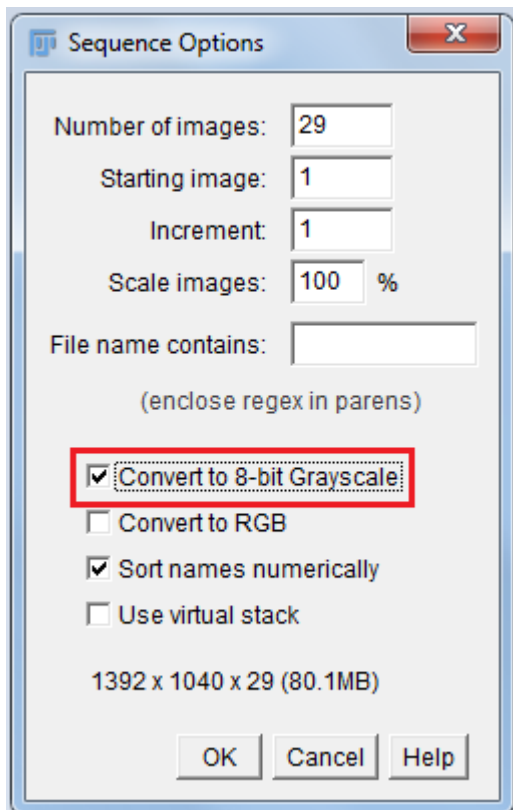


- Lunch single-cell analysis plugin *Plugins->MASTER PLUGIN*
- Select following options before running tool

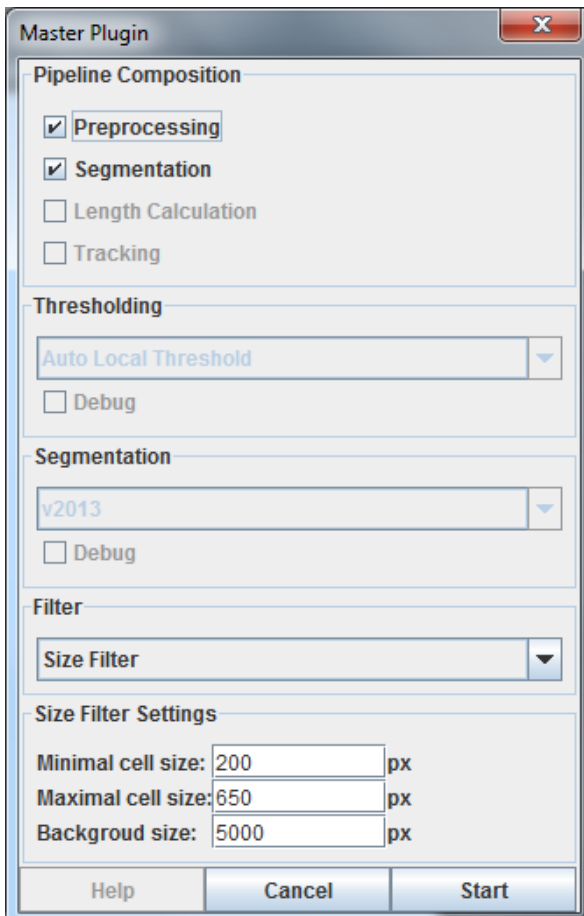

- The tool is trying to automatically detect the cultivation chamber, if it couldn't be detected you have to select it manually.

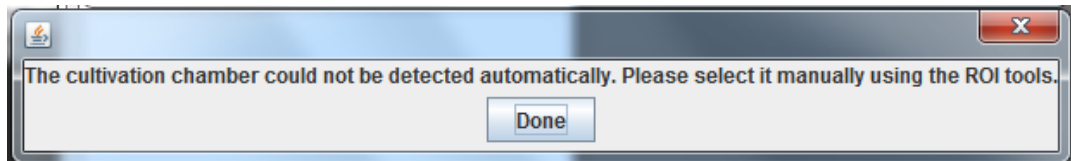

- Use the rectangular selection tool.

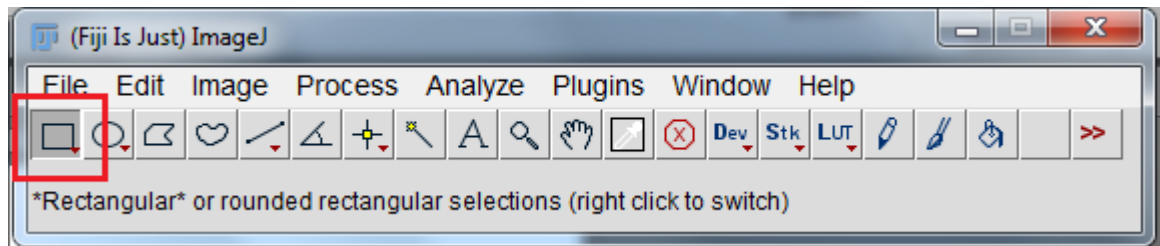

- To select the cultivation chamber as such:




- Afterwards click done and the plugin is going to detect and analysis the single cells in the frames

Post-processing

After a couple of seconds the tool has detected various objects, not always cells, in each of the given frames. The detected objects are highlighted by yellowish curves and are listed in the *ROI (region of interest) Manager*. Furthermore, the number of cells and overall size of each frame is shown in the *General Results* windows.



Imported images series

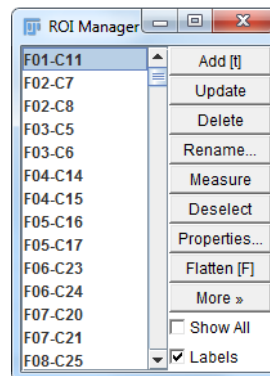

File Edit Font		
Frame	Cell number	Overall size
1.000	1.000	637.000
2.000	2.000	715.000
3.000	2.000	800.000
4.000	2.000	891.000
5.000	2.000	1003.000
6.000	2.000	1110.000
7.000	2.000	1197.000
8.000	3.000	1309.000
9.000	4.000	1486.000
10.000	4.000	1668.000


ROI Manager and General Results windows

Go through each of the frames and check if:

- only cells are detected and no other objects




Wrong detected object

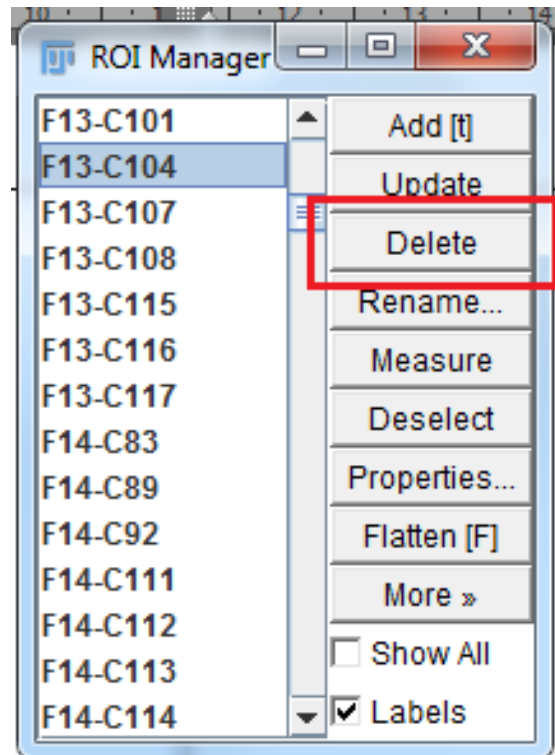

Delete wrong detected object

- that the cell was recognized entirely

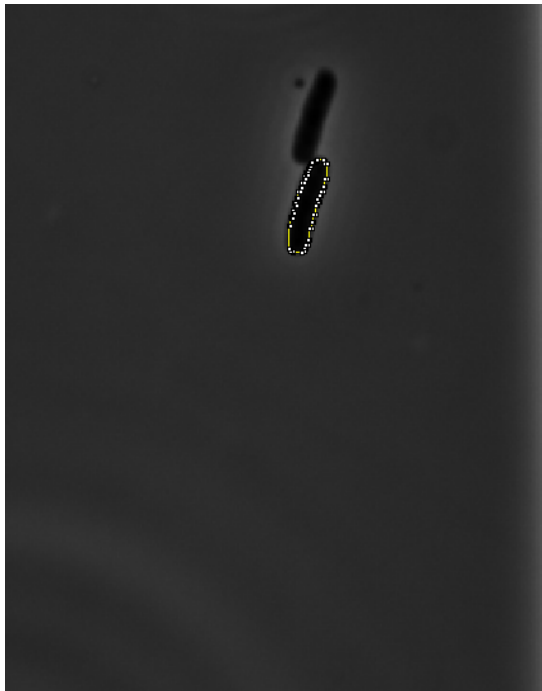


Cell was detected entirely, yellowish curve surround it completely.




Cell was not detected entirely, yellowish curve doesn't surrounds it completely.

Use the *Selection Brush Tool* for altering the ROI, use the shortcut *shift + left mouse click* to enlarge the ROI to fit the entire cell shape. Merge ROI by pressing **m** on the keyboard.

By double clicking on the *Selection Brush Tool* let you change the size of the brush in pixels.

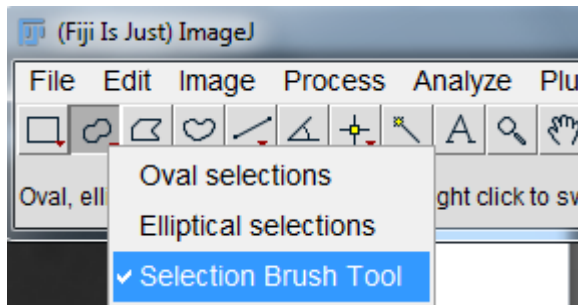


Selection of brush tool for manual alteration of ROI area

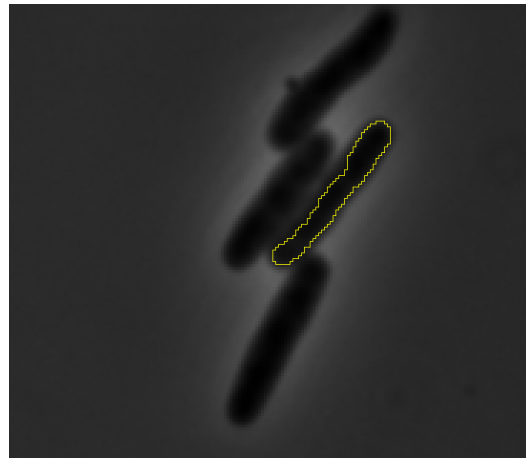


After altering the ROI area with the brush tool press m to make the change permanent

- cells were separated correctly:


Use again the *Selection Brush Tool* and the key combination *control + left mouse click* to split the ROI in two separate ones. Split ROIs by pressing **s** on the keyboard.

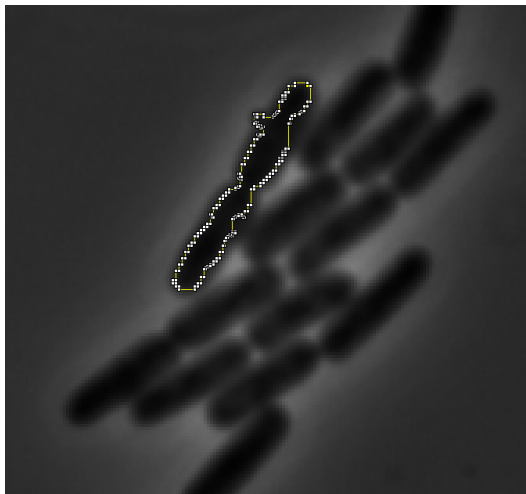


Not proper separated cells

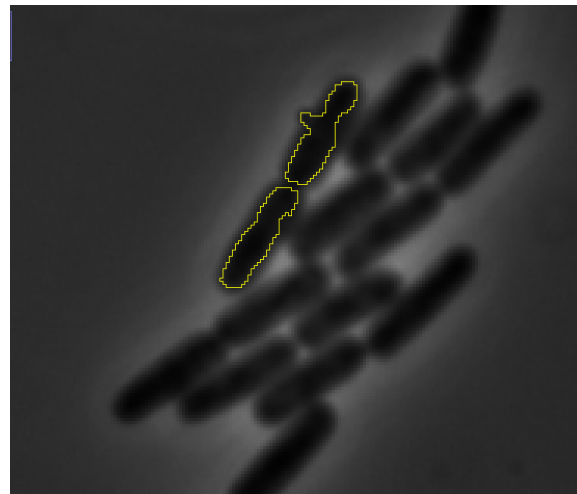

Separated ROIs using the *Selection Brush Tool*. Press s on the keyboard to make the change permanent.

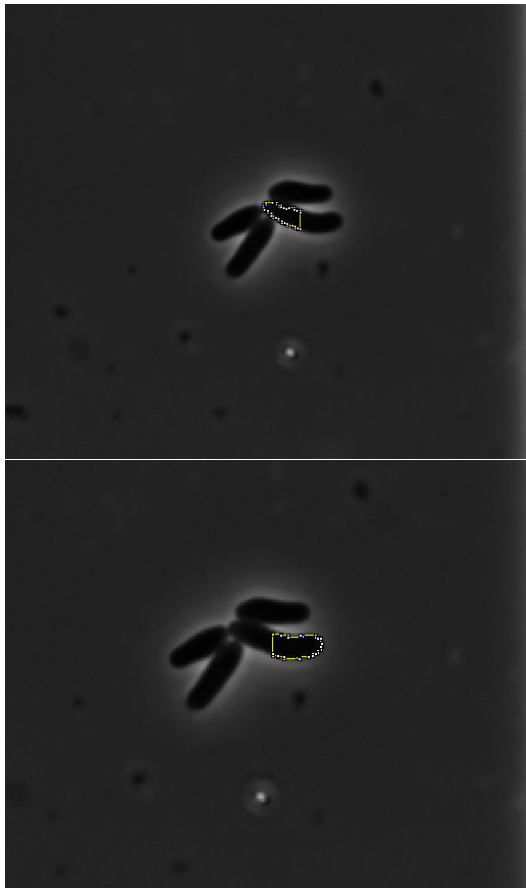

Cell was wrongly separated in two

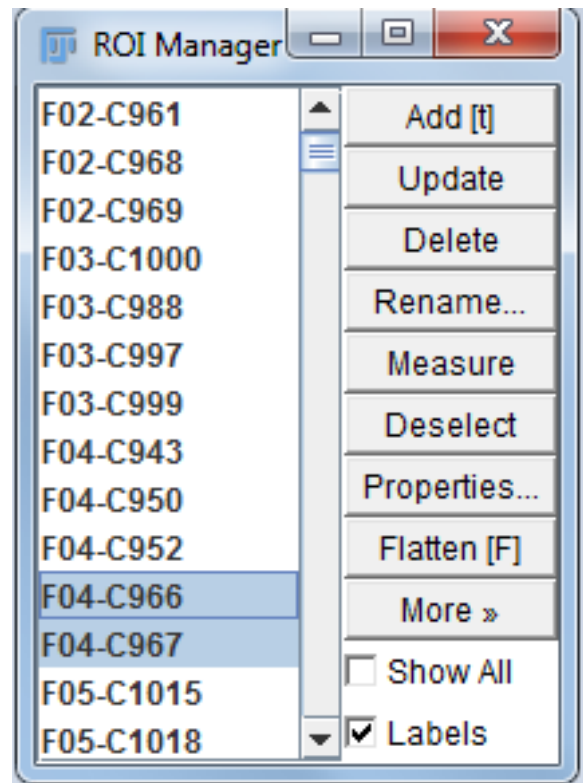

By selecting both ROIs in the ROI Manager both can be rejoined by pressing **m** on the keyboard


- Select both ROIs in the RIO Manager and press **m** to merge them

Measurement data

Save analysed images and ROIs first

- Save the tif image series as zip file using *File->Save as->ZIP...*
- Select all ROIs in the ROI Manger using the key combination *ctrl + a* and save them *More -> Save ...*

Measure ROIs and export results

- Click on  to create a new *General Results* table.
- Save data using *File->Save as ...*

Exemplary data sets

You can find exemplary data sets containing pre-processed and segmented image series of *Corynebacterium glutamicum* and *Escherichia coli* here:

Troubleshooting

Please report any issues regarding the plugin here: <http://ibtmodesimhub/trac/ij-plugin/>

You can log in with your normal user and issue a new ticket. Please try to describe your problem as detailed as possible.